

**ICONIX™**  
All Suture Anchor

with IntelliBraid Technology




# Introduction

## Searching For a Better Option


The ICONIX family of anchors symbolize the next generation of suture anchor technology. The all suture based system allows for less bone removal during pilot hole creation without compromising pullout strength. The innovative design is truly the first of its kind.

## IntelliBraid Technology

The IntelliBraid Technology creates a bunching effect using targeted compression zones within the implant sheath for optimal fixation.

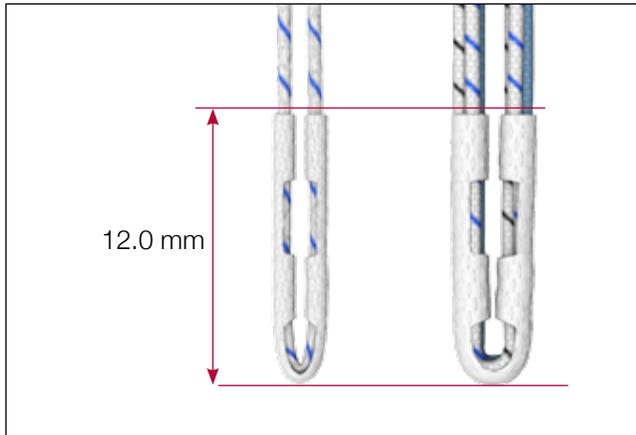

Measurements - Before

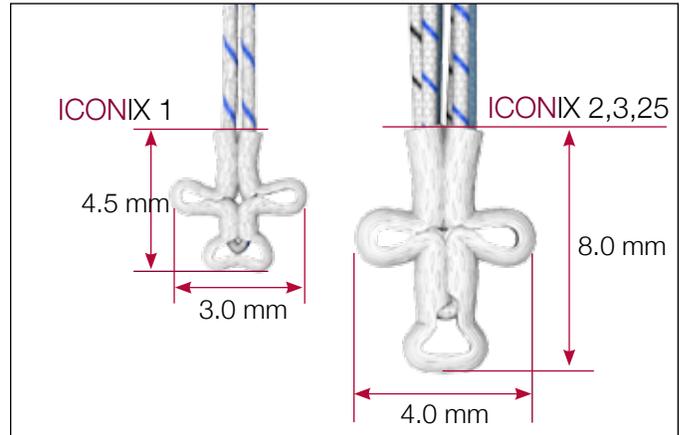

Measurements - After

## Self-Centering Technology

The ICONIX disposable drills feature a unique Self Centering (SC) Technology, to ensure accurate pilot hole placement.



Without SC Technology



With SC Technology

# Competitive Advantages

## Anchors

- Four all suture configurations
- A completely all suture based implant made of HSS (Force Fiber) and a polyester sheath
- Helps eliminate the possibility of loose bodies in the joint due to anchor breakage



**ICONIX 1**  
1 strand #2 Force Fiber  
Use w/1.4mm instrumentation



**ICONIX 2**  
2 strands #2 Force Fiber  
Use w/2.3mm instrumentation




**ICONIX 3**  
3 strands #2 Force Fiber  
Use w/2.3mm instrumentation



**ICONIX 25**  
2 strands #5 Force Fiber  
Use w/2.3mm instrumentation

## Pilot Hole

- Smaller drill allows for minimal bone removal during pilot hole creation
- More anchors can be placed in smaller footprint

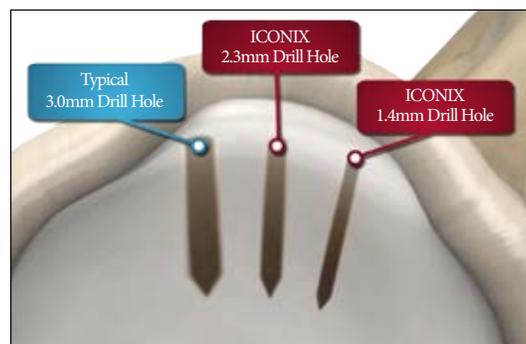

## Intelligent Design

- The ICONIX Anchor System is intelligently sized and color coded to reduce risk of human error

**1.4mm**



**2.3mm**



# Competitive Advantages

## Guides

- Smaller drill guides cause less damage to surrounding tissue
- Universal crown tips can be used on the face and rim of the glenoid



- Options when you need them
- The ICONIX Anchor System was specifically designed for both Rotator Cuff and Instability repairs, with curved and straight option guides

### *Rotator Cuff*

— — — — — Straight Short — — — — —

**2.3mm**


### *Instability*

Straight

**1.4mm**


12°

25°

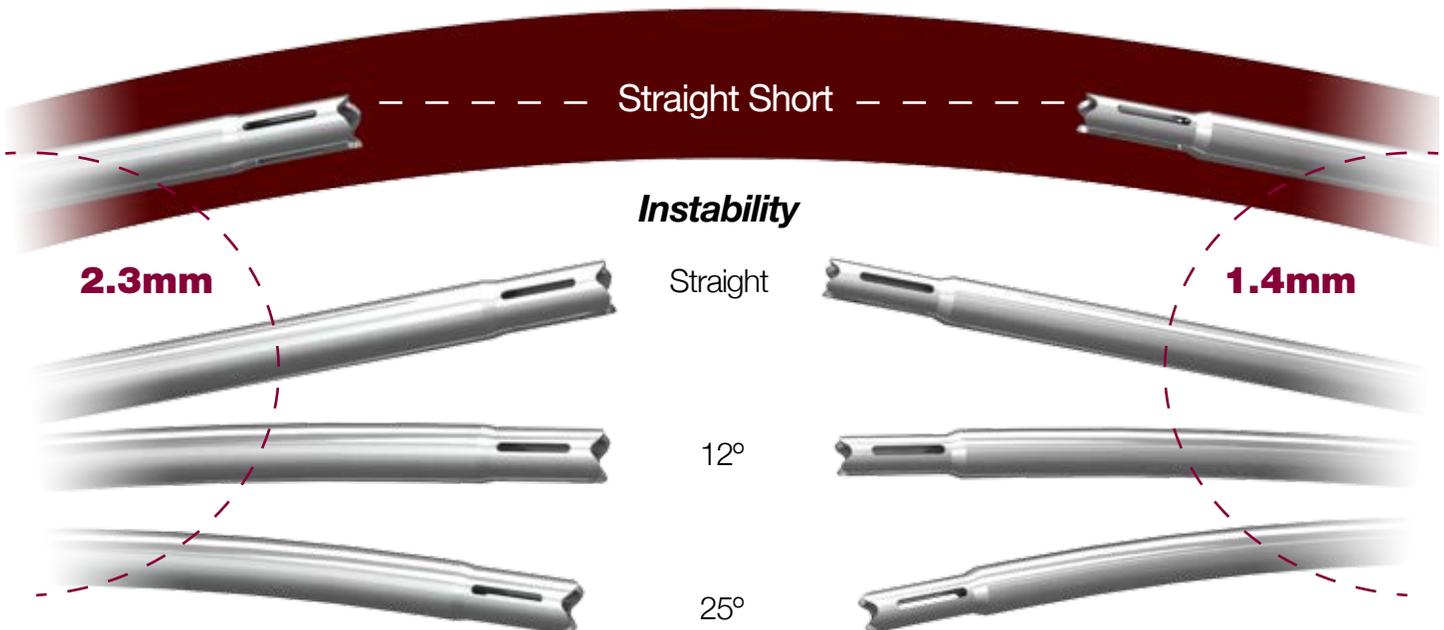


# Competitive Advantages

- The ICONIX Anchor System is specifically designed for Hip FAI procedures with curved and straight guide options
- Designed to provide proper guide length during all FAI cases

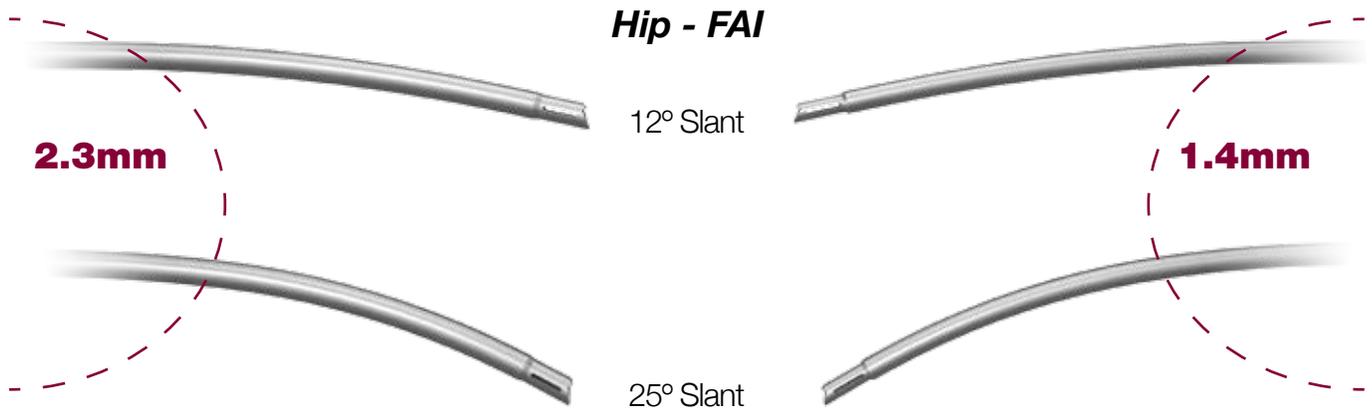

- Directs anchor trajectory away from the acetabular socket
- Unique SLANT mouth helps place anchor closer to the acetabulum rim, allowing for a more anatomic repair

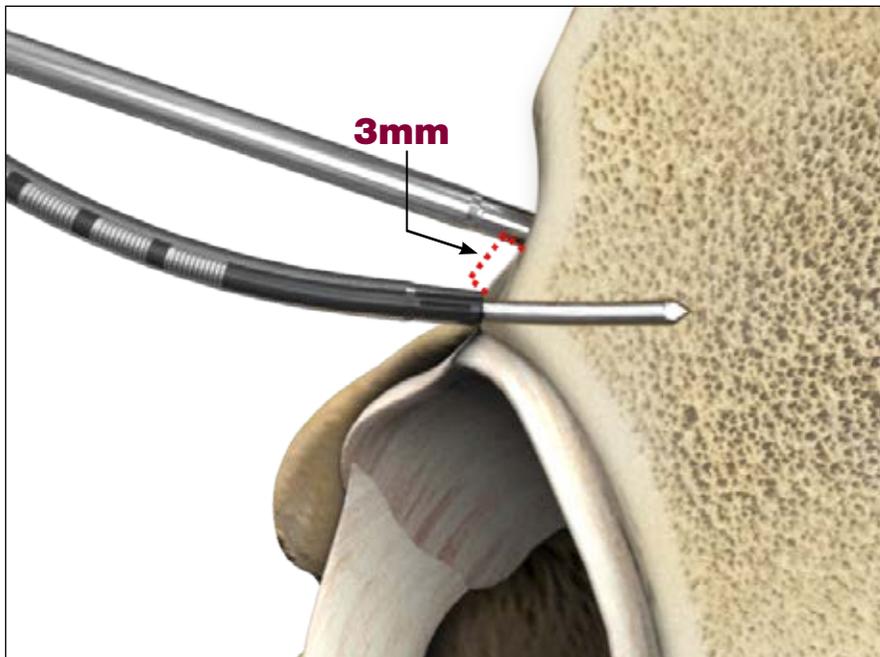


# Technique


## Rotator Cuff



Take selected drill guide and place on bone in preferred position.

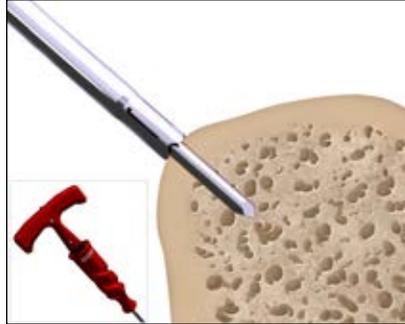

Insert the awl or drill into the drill guide and advance until positive stop is hit.

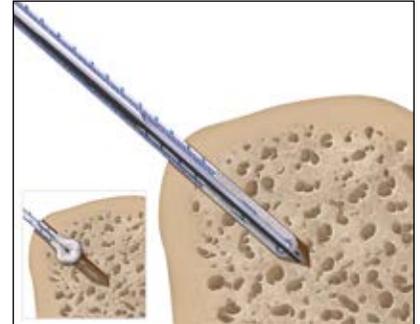


Insert anchor, lightly mallet to fully seat the anchor into the bone until positive stop is hit. Untether the suture from the inserter handle. Pull the anchor inserter handle. Pull on all sutures to set the anchor.


## Instability



Take selected drill guide and place on bone in preferred position.

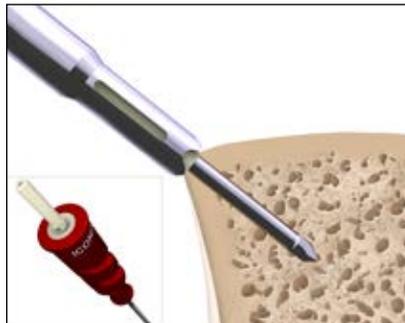

Insert the drill into the drill guide and advance until positive stop is hit.

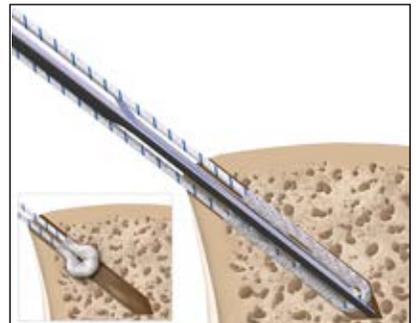


Insert anchor, lightly mallet to fully seat the anchor into the bone until positive stop is hit. Untether the suture from the inserter handle. Pull the anchor inserter handle. Pull on all sutures to set the anchor.

## Hip - FAI



(Instructions for Hip - FAI are same as for Instability)

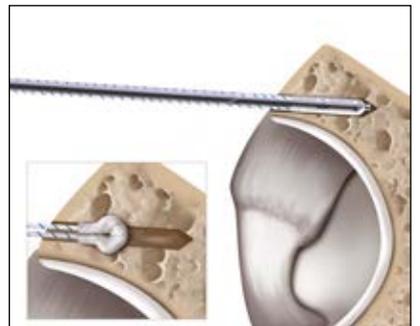


# Competitive Advantages

## Cannulated Obturator

- Cannulated obturators help save time while gaining access to the hip joint



## Pull-out Strength<sup>1</sup>

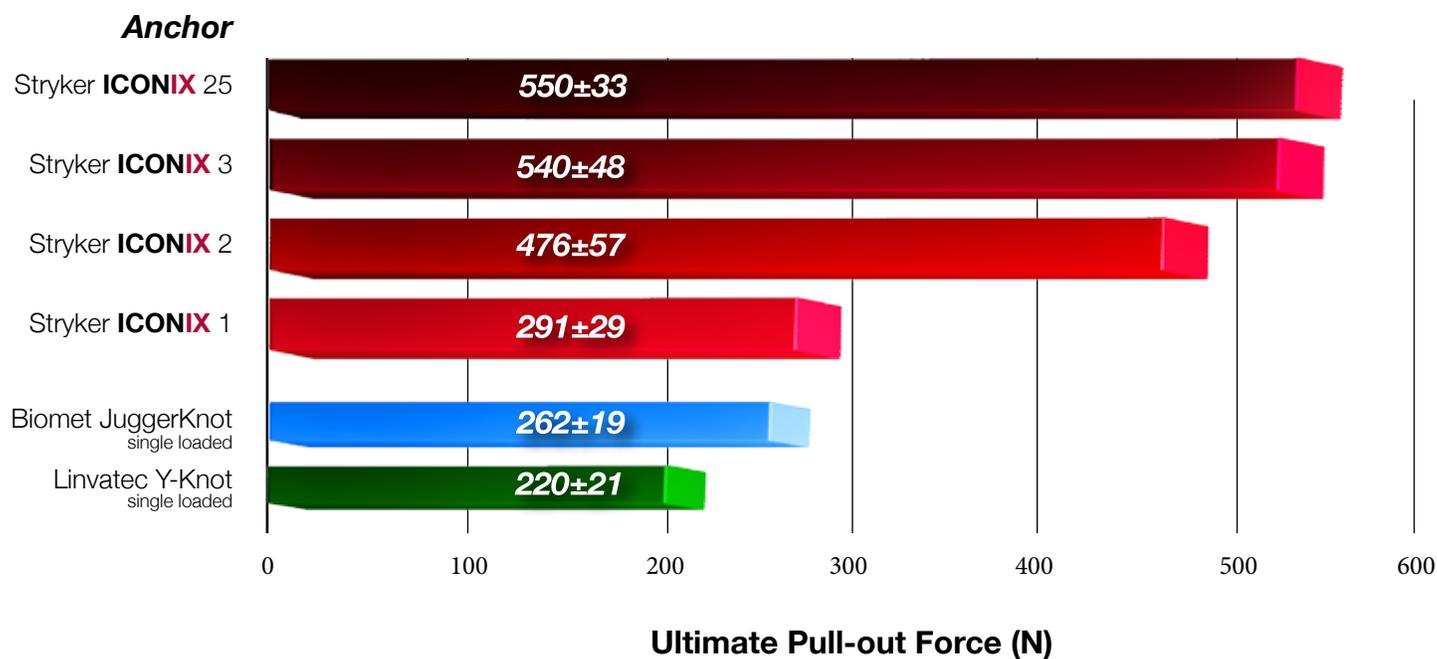

One pound force is approximately equal to 4.5 Newtons

<b>PART NUMBER</b>	<b>DESCRIPTION</b>
<b>Anchors</b>	
3910-500-512	ICONIX 1 with IntelliBraid Technology, 1.4mm Anchor with 1 strand #2 Force Fiber
3910-500-522	ICONIX 2 with IntelliBraid Technology, 2.3mm Anchor with 2 strands #2 Force Fiber
3910-500-532	ICONIX 3 with IntelliBraid Technology, 2.3mm Anchor with 3 strands #2 Force Fiber
3910-500-525	ICONIX 25 with IntelliBraid Technology, 2.3mm Anchor with 2 strands #5 Force Fiber
<b>Guides</b>	
3910-500-550	ICONIX Straight Guide SHORT for 1.4mm anchor
3910-500-551	ICONIX Straight Guide SHORT for 2.3mm anchor
3910-500-552	ICONIX Straight Guide for 1.4mm anchor
3910-500-553	ICONIX 12' Guide for 1.4mm anchor
3910-500-554	ICONIX 25' Guide for 1.4mm anchor
3910-500-555	ICONIX Straight Guide for 2.3mm anchor
3910-500-556	ICONIX 12' Guide for 2.3mm anchor
3910-500-557	ICONIX 25' Guide for 2.3mm anchor
3910-500-558	ICONIX SLANT 12' Guide for 1.4mm anchor
3910-500-559	ICONIX SLANT 25' Guide for 1.4mm anchor
3910-500-560	ICONIX SLANT 12' Guide for 2.3mm anchor
3910-500-561	ICONIX SLANT 25' Guide for 2.3mm anchor
<b>Obturator</b>	
3910-500-562	ICONIX Pencil Tip Obturator SHORT for 1.4mm anchor
3910-500-563	ICONIX Pencil Tip Obturator SHORT for 2.3mm anchor
3910-500-564	ICONIX Trocar Tip Obturator for 1.4mm anchor
3910-500-565	ICONIX Trocar Tip Obturator for 2.3mm anchor
3910-500-566	ICONIX Cannulated Bullet Obturator for 1.4mm anchor
3910-500-567	ICONIX Cannulated Bullet Obturator for 2.3mm anchor
<b>Drills / Awl</b>	
3910-500-568	ICONIX 1.4mm Disposable Drill
3910-500-569	ICONIX 2.3mm Disposable Drill
3910-500-571	ICONIX 1.4mm Reusable Straight Drill
3910-500-573	ICONIX 2.3mm Reusable Straight Drill
3910-500-574	ICONIX 2.3mm Reusable Awl
<b>Tray</b>	
3910-500-575	Multi-System Reusable Tray

References:  
Tech Report # RD-12-080

A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any Stryker product. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: ICONIX, Stryker. All other trademarks are trademarks of their respective owners or holders.

Force Fiber is a registered trademark of Teleflex Incorporated.

Literature Number: LJPICX-B  
MS/GS 09/12

Copyright © 2012 Stryker  
Printed in USA

325 Corporate Drive  
Mahwah, NJ 07430  
t: 201 831 5000

[www.stryker.com](http://www.stryker.com)