



Versatility & Control, Performance & Reliability

## InterPulse Battery Powered Irrigation

### InterPulse Handpiece

PART #	DESCRIPTION
210-100	Interpulse Handpiece with suction
210-110	Interpulse Handpiece with suction with Bone Cleaning Tip
210-110-100	Interpulse Handpiece with suction with Coaxial Bone Cleaning Tip
210-112	Interpulse Handpiece with suction with Soft Tissue Tip
210-114	Interpulse Handpiece with suction with High Flow Tip
210-114-100	Interpulse Handpiece with suction with Coaxial High Flow Tip
210-118	Interpulse Handpiece with suction with Fan Spray Tip
210-118-100	Interpulse Handpiece with suction with Coaxial Fan Spray Tip
210-200	Irrigation only Interpulse Handpiece
210-210	Irrigation only Interpulse Handpiece with Bone Cleaning Tip
210-210-100	Irrigation only Interpulse Handpiece with Coaxial Bone Cleaning Tip
210-212	Irrigation only Interpulse Handpiece with Soft Tissue Tip
210-214	Irrigation only Interpulse Handpiece with High Flow Tip
210-214-100	Irrigation only Interpulse Handpiece with Coaxial High Flow Tip
210-218	Irrigation only Interpulse Handpiece with Fan Spray Tip
210-218-100	Irrigation only Interpulse Handpiece with Coaxial Fan Spray Tip

### Tips and Accessories

PART #	DESCRIPTION
210-1	Wound Tunneling Tip with Suction
210-3	Wound Tunneling Tip
210-4	Femoral Canal Brush
210-5	Femoral Canal Tip (Irrigation Only)
210-7	Coaxial Femoral Canal Tip
210-8	Femoral Canal Tip
210-10	Bone Cleaning Tip
210-10-100	Coaxial Bone Cleaning Tip
210-12	Soft Tissue Tip
210-14	High Flow Tip
210-14-100	Coaxial High Flow Tip
210-16	Irrigation Only Short Tip
210-17	Revision Tip
210-18	Fan Spray (Low Splash) Tip
210-18-100	Coaxial Fan Spray Tip
210-34	Medium Soft Cone Splash Shield
210-312	Soft Tissue Tip with Replacement Suction Tubing
210-318	Fan Spray Tip with Replacement Suction Tubing
211-22	Interpulse Rechargeable Power Pack
211-24	Interpulse Bay Charger
211-100	Interpulse Rechargeable Hand Piece

### InterPulse Handpiece, Rechargeable

PART #	DESCRIPTION
211-100	Interpulse Rechargeable Handpiece
211-110	Interpulse Rechargeable Handpiece with Bone Cleaning Tip
211-110-100	Interpulse Rechargeable Handpiece with Coaxial Bone Cleaning Tip
211-112	Interpulse Rechargeable Handpiece with Soft Tissue Tip
211-114	Interpulse Rechargeable Handpiece with High Flow Tip
211-114-100	Interpulse Rechargeable Handpiece with Coaxial High Flow Tip
211-118	Interpulse Rechargeable Handpiece with Fan Spray Tip
211-118-100	Interpulse Rechargeable Handpiece with Coaxial Fan Spray Tip
211-200	Interpulse Rechargeable Irrigation Only Handpiece
211-210	Interpulse Rechargeable Irrigation Only Handpiece with Bone Cleaning Tip
211-210-100	Interpulse Rechargeable Irrigation Only Handpiece with Coaxial Bone Cleaning Tip
211-212	Interpulse Rechargeable Irrigation Only Handpiece with Soft Tissue Tip
211-214	Interpulse Rechargeable Irrigation Only Handpiece with High Flow Tip
211-214-100	Interpulse Rechargeable Irrigation Only Handpiece with Coaxial High Flow Tip
211-218	Interpulse Rechargeable Irrigation Only Handpiece with Fan Spray Tip
211-218-100	Interpulse Rechargeable Irrigation Only Handpiece with Coaxial Fan Spray Tip

\* Diane Morgan, BSN, RNC, CWCN, CETN  
 Director of Clinical Services  
 LifeSpan Home & Hospice Care, Battle Creek, MI



For more information, contact your local Stryker Sales Representative.

4100 East Milham Ave  
 Kalamazoo, MI 49001  
 (800)253-3210  
 (616)323-7700

**stryker**  
**INSTRUMENTS**

# InterPulse

Battery Powered Irrigation



# Wound Care

# Performance, Reliability, Versatility

Use Interpulse to prepare the wound bed to heal the wound more quickly.

## Features/Benefits

- Ergonomic shape to reduce muscle fatigue
- FDA approved PSI <15
- Single patient use
- Portable
- Cost effective



**Tip Locking Mechanism**  
Effectively engages tip to hand piece for optimal performance.

**Full Speed Lock Mechanism**  
The Interpulse can lock in at full speed.

**Simple to Use**  
Simply spike the battery pack into the saline bag, insert tip into hand piece and pull the trigger.

**Lightweight / Portable**  
Battery pack holds 8 AA batteries attached to the spike for light weight handpiece.

## Effective Wound Debridement



Chronic venous insufficiency before hydrodebridement treatment with InterPulse.



Wound free of necrotic tissue after 5 treatments. Follow up treatments continue to healthy granulating wound.

Photos courtesy of Diane Morgan  
Llispian Home & Hospice Care



# Wound Care Replacement Tubing: Step-by-Step


Reusing your InterPulse handpiece with the Stryker replacement tubing kits



## Step-by-Step



2. Remove old tubing by pulling up and out of top of handpiece.  
3. Pinch new tubing and press into channel on top rear of the handpiece.  
4. Pull tubing forward and down into channel while gripping handpiece firmly.  
5. Pull tubing to the rear or cut so tubing is flush with distal end of handpiece.  
6. Tubing should be flush with distal end of handpiece for proper tip insertion.

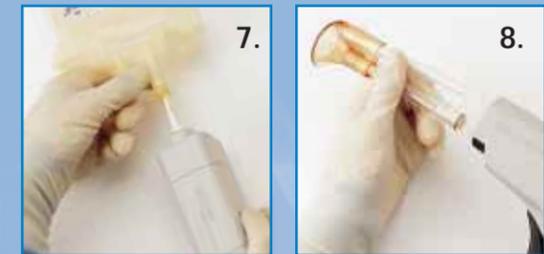

7. Spike saline bag.  
8. Insert tip until it "clicks".

210-3 Wound Tunneling Tip  
210-1 Wound Tunneling Tip with Suction



210-18-100 Coaxial Fan Spray Tip with Pliable Splash Shield  
210-18 Fan Spray Tip  
210-12 Soft Tissue Tip

Katie, would it be appropriate to fill this space with a brief note on the flexible wound tunneling tips? Just a suggestion.

**Advanced tip technology & selection for your Wound Care needs.**