



Signature Portfolio MIS Attachments and Cutting Accessories

Attachments	
Part Number	Product Description
5407-120-900	MIS Hub
5407-120-950	13cm MIS Straight Nose Tube Color Band: Lilac
5407-120-950A	13cm MIS Angled Nose Tube Color Band: Lilac
5407-120-950C	13cm MIS Curved Nose Tube Color Band: Lilac
5407-120-970	16cm MIS Straight Nose Tube Color Band: Brown
5407-120-970A	16cm MIS Angled Nose Tube Color Band: Brown
5407-120-970C	16cm MIS Curved Nose Tube Color Band: Brown

To order MIS attachments and cutting accessories – or to learn more about our Signature Portfolio – please visit www.stryker.com/NSE or call **800-253-3210**.

NSE Neuro Spine ENT

We're laser focused on Neuro, Spine and ENT. Our instruments, technology, education and support are all tailored to your specialties to help fulfill your unique needs and potentially make a positive difference in your day and outcomes.

- High Speed Drills (electric and pneumatic)
- Attachments and cutting accessories
- Sonopet® Ultrasonic Aspirator
- Bone Mill
- SilverGlide® Bipolar Forceps
- Microdebrider
- NasoPore® Bioresorbable Nasal Dressings
- HemoPore® Hemostatic Bioresorbable Nasal Dressing
- OtoPore® Bioresorbable Ear Dressings
- Surgeon headlights and other accessories in support of NSE specialties

Cutting Accessories		
Part Number	Product Description	Color Band
8450 (lilac) series = 13cm shank 8470 (brown) series = 16cm shank		
8450-107-015	1.5mm Match Head	Lilac
8470-107-015	1.5mm Match Head	Brown
8450-107-020	2.0mm Match Head	Lilac
8470-107-020	2.0mm Match Head	Brown
8450-107-025	2.5mm Match Head	Lilac
8470-107-025	2.5mm Match Head	Brown
8450-107-030	3.0mm Match Head	Lilac
8470-107-030	3.0mm Match Head	Brown
8450-107-525	2.5mm Precision Match Head	Lilac
8470-107-525	2.5mm Precision Match Head	Brown
8450-107-530	3.0mm Precision Match Head	Lilac
8470-107-530	3.0mm Precision Match Head	Brown
8450-107-015D	1.5mm Diamond Match Head	Lilac
8470-107-015D	1.5mm Diamond Match Head	Brown
8450-107-020D	2.0mm Diamond Match Head	Lilac
8470-107-020D	2.0mm Diamond Match Head	Brown
8450-107-025D	2.5mm Diamond Match Head	Lilac
8470-107-025D	2.5mm Diamond Match Head	Brown
8450-107-030D	3.0mm Diamond Match Head	Lilac
8470-107-030D	3.0mm Diamond Match Head	Brown
8450-107-035D	3.5mm Diamond Match Head	Lilac
8470-107-035D	3.5mm Diamond Match Head	Brown
8450-107-040D	4.0mm Diamond Match Head	Lilac
8470-107-040D	4.0mm Diamond Match Head	Brown
8450-107-030DC	3.0mm Coarse Diamond Match Head	Lilac
8470-107-030DC	3.0mm Coarse Diamond Match Head	Brown
8450-107-040DC	4.0mm Coarse Diamond Match Head	Lilac
8470-107-040DC	4.0mm Coarse Diamond Match Head	Brown
8450-010-020	2.0mm Round	Lilac
8470-010-020	2.0mm Round	Brown
8450-010-030	3.0mm Round	Lilac
8470-010-030	3.0mm Round	Brown
8450-010-040	4.0mm Round	Lilac
8470-010-040	4.0mm Round	Brown
8450-009-020	2.0mm Precision Round	Lilac
8470-009-020	2.0mm Precision Round	Brown
8450-009-030	3.0mm Precision Round	Lilac
8470-009-030	3.0mm Precision Round	Brown
8450-009-040	4.0mm Precision Round	Lilac
8470-009-040	4.0mm Precision Round	Brown
8450-009-050	5.0mm Precision Round	Lilac
8470-009-050	5.0mm Precision Round	Brown
8450-012-020DC	2.0mm Coarse Diamond	Lilac
8470-012-020DC	2.0mm Coarse Diamond	Brown
8450-012-030DC	3.0mm Coarse Diamond	Lilac
8470-012-030DC	3.0mm Coarse Diamond	Brown
8450-012-040DC	4.0mm Coarse Diamond	Lilac
8470-012-040DC	4.0mm Coarse Diamond	Brown
8450-013-040DX	4.0mm Diamond Extra Coarse	Lilac
8470-013-040DX	4.0mm Diamond Extra Coarse	Brown
8450-013-050DX	5.0mm Diamond Extra Coarse	Lilac
8470-013-050DX	5.0mm Diamond Extra Coarse	Brown

stryker®

Reconstructive

Hips

Knees

Trauma & Extremities

Joint Preservation

Orthobiologics

Medical & Surgical

Power Tools & Surgical Accessories

Image Guided Navigation

Endoscopy & Arthroscopy

Integrated Communications

Beds, Stretchers & EMS

Sustainability Solutions

Neurotechnology & Spine

Craniomaxillofacial

Interventional Spine

Neurosurgical, Spine & ENT

Neurovascular

Spinal Implants

The information presented in this brochure is intended to demonstrate the breadth of Stryker product offerings. Always refer to the package insert, product label and/or user instructions before using any Stryker product. Products may not be available in all markets. Product availability is subject to the regulatory or medical practices that govern individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Products referenced with TM designation are trademarks of Stryker. Products referenced with ® designation are registered trademarks of Stryker.

Literature Number: 9100-003-174 Rev None G55/PS

Copyright © 2015 Stryker
Printed in USA

Stryker Neuro Spine ENT
4100 East Milham Avenue
Kalamazoo, MI 49001 USA
t: 269 323 7700 f: 800 999 3811
toll free: 800 253 3210

www.stryker.com/NSE

stryker®

Neuro Spine ENT

Signature
PORTFOLIO

MIS Attachments



MIS

Maximizing functionality and options, while being mindful of cost

Our Signature Portfolio high speed drill platform offers a robust, new array of minimally invasive surgery (MIS) attachments and cutting accessories. Based on user input, attachments have been rebalanced as well as refined in multiple ways to enhance your hands-on experience. Cutting accessories have also been optimized, giving you the largest selection of MIS cutting accessories and the largest MIS cutting accessory head on the market today.

Each MIS cutting accessory provides multiple choices for the precise reach, customization and exposure you want for a particular patient or procedure. Large diameter cutting accessories offer four different telescoping exposure points (in 1mm increments) and small diameter cutting accessories offer ten; all easily adjusted using built-in notches for exposure and minimum/maximum lengths. This flexibility to extend and retract each cutting tool allows you to cost effectively secure the broad options and functionality you want through a streamlined set of instruments in your inventory or surgical field.

Here's a closer look at our two series of MIS attachments (13 and 16cm), and their corresponding cutting accessories.



The only MIS cutting accessories on the market with multi-notch (4-10mm) telescoping ability; maximizes length options while minimizing cost and inventory

Refined attachment design for improved balance and maneuverability

Alignment dots and Snap-Connect guide assembly/disassembly of attachments to motor in an easy and consistent manner

Color coding instantly identifies matching attachments and cutting accessories

Longer, concave design for improved line of sight and fine-edge knurling for better grip and more options for positioning your hand

▲ Patented 2Flute cutting accessory, smooth and precise for the rapid dissection of bone

▲ Head diameter up to 5mm, largest on the market

Portfolio options

Signature Portfolio welcomes more than 50 new MIS cutting accessories, examples of which are shown below. To suit your every need, they come in an unparalleled variety of head types, diameters (1.5 - 5.0mm) and lengths (13-16cm). Please see the back cover for our full lineup.

Match Head



Precision Match Head




Diamond Match Head




Coarse Diamond Match Head




Round

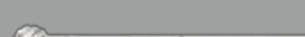

Precision Round

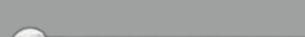

Coarse Diamond



Diamond Extra Coarse

